



BioFactory™

CREATE THREE-DIMENSIONAL ORGANOMIMETIC MODELS FOR TISSUE ENGINEERING








BIOPRINTING TECHNOLOGY IS EXPANDING IN RESEARCH AND INDUSTRY

BioFactory™

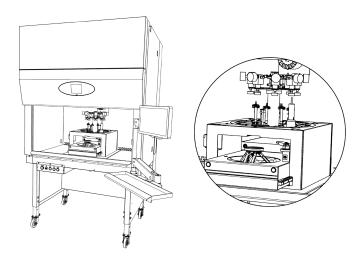

> The BioFactory™ is a versatile and cell friendly three-dimensional manufacturing instrument allowing scientists to pattern cells, biomolecules and a range of soft and rigid materials in desirable 3D composite structures, mimicking natural environments.

GENERAL SPECIFICATIONS	
Working range	$60 \times 60 \times 60$ mm
High degree of customization	
Precision	± 5 μm
Modular printhead concept	GENERAL
Nano liter dispensing resolution, minir	mal dead volume
Printing under physiological conditions	s
Temperature control	from 5°C up to 80°C
	(substrate holder, medias)
Overall dimensions (W $ imes$ L $ imes$ H)	1370 imes 1030 imes 2400 mm

A POWERFUL TISSUE ENGINEERING INSTRUMENT TO MIMIC NATURE

> It has been well documented that cell culture in tissue culture plates and flasks doesn't mimic the in-vivo cell growth and has poor correlation with in-vivo animal as well as human clinical data.

The BioFactory™ instrument provides a powerful tool for tissue engineering to create organotypic tissues with in-vivo-like morphology that better mirrors the environment experienced by cells in-vivo. Furthermore, it better reflects cell behavior, intercellular interactions and differentiation processes.

BIOPRINTING PLATFORM

- > The BioFactory™ instrument is a cutting-edge platform to explore the potential of 3D tissue engineering through bioprinting technology.
- Spatial control of cells, bioactives, and extracellular matrix in a threedimensional cellular construct is an enabling approach to construct designed organomimetic tissues for drug discovery and regenerative medicine.

KEY ADVANTAGES

- $> \ \ \, \text{Printing under physiological conditions}$
- > Laser units for photo polymerization/ biomolecule immobilization
- > Biomanufacturing within a sterile laminar flow hood
- > Repeatability of micrometer-scale processes
- $>\;\;$ Fast and easy tissue modelling via Bioprinting Suite
- > Customization options



TOOLS, OPTIONS AND ACCESSORIES

Laser or photo-crosslinking devices:	for hydrogel polymerization, bioactives encapsulation, signal molecule immobilization,	
	coating or ablation processes	
Printhead technologies:	for optimal processing of abroad biomaterial/bioactives portfolio:	
	 cell-friendly lnk-jet; 	
	- thermopolymer extrusion;	
	 2 component printhead; 	
	 paste and hydrogel dispensing. 	
High precision temperature control devices:	for biomaterials, mediums, printheads and substrates	
Electrospinning printhead:	printhead for micro & submicrometer bioarchitectures manufacturing	
Software Suite:	to interact with bioprinting instrument, medical imaging, human machine interface	
	 including BioCAD™, BioCAM™, BioCUT™, Industry standard interfaces (STL, DICOM, AMF, DXF) 	
Calibration systems:	calibration laser, needle, substrate / lab device	



SWISS INNOVATION